



Forma® ADAPT Splint

Hard function, soft feeling



Trees for All



Elysee Dental (NL) | Flemingweg 14 | 2408 AV Alphen aan den Rijn |
T +31 (0)172 65 33 00 | info@elysee-dental.nl | www.elysee-dental.nl



Hard function, soft feeling

A splint does not treat the clenching or grinding activity, but it protects the teeth and jaw system against overload.


The clinical data from your IOS records are used to design the desired occlusal device digitally, then 3D printed using **Keysplint Soft™** material.

This revolutionary material is hard enough for a functional occlusal device, while soft enough to provide your patient with the desired comfort and durability.

** Keysplint Soft™ is a brand by **keystone***

Any design is possible

- Stabilisation (Michigan, without clamps).
- Disclusion hills (canine region), 1 mm lateral freedom of shift.
- Concave frontal velodrome.
- Protection.
- Occlusal thickness of 1,5 - 2 mm.
- Peripheral thickness of 1,5 - 2 mm.
- Flat occlusal (no imprints).

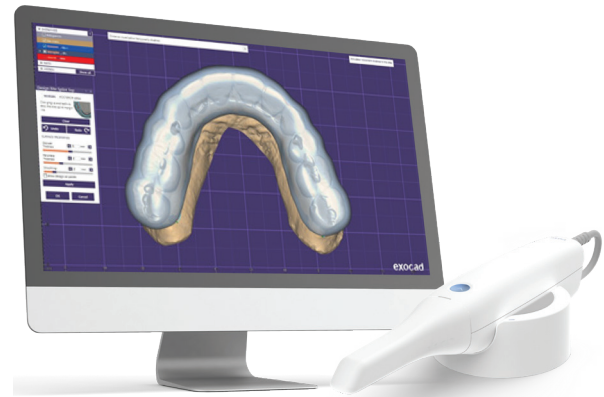

Only for IOS cases

For the manufacturing of the Forma® ADAPT splint, only three scans are required:

- Maxillary arch
- Mandibular arch
- VDO bite registration (leaf gauges technique)

Benefits of the Forma® ADAPT

- Digital accuracy.
- Strength & flexibility.
- Memory function.
- Ability to reprint from design.
- Biocompatible material.
- Enhanced patient comfort.
- Stain & fracture resistant.
- No model necessary.

